

ENGLISH

For Professional Use Only

KMT2A-MLLT1 Dual Fusion/Translocation FISH Probe Kit

Introduction


The KMT2A-MLLT1 Fusion/Translocation FISH Probe Kit is designed to detect rearrangements involving the human KMT2A and MLLT1 genes located on chromosome bands 11q23.3 and 19p13.3, respectively. Rearrangements between the two genes, the KMT2A gene - also known as HRX, MLL, MLL1, TRX1, ALL-1, CXXC7, HTRX1, MLL1A, WDSTS, MLL/GAS7 or TET1-MLL – and the MLLT1 gene – also called ENL, LTG19 or YEATS1, have been observed in acute myeloid leukemia (AML), acute lymphoblastic leukemia (ALL) and other malignancies.

Intended Use

To detect rearrangements involving the human *KMT2A* and *MLLT1* genes located on chromosome bands 11q23.3 and 19p13.3, respectively.

Cont.	Color
LSP KMT2A 5'-3' FISH Probe	CytoOrange
LSP MLLT1 5'-3' FISH Probe	CytoGreen

Probe Design

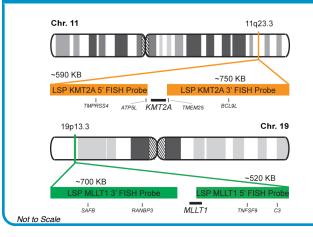

LSP KMT2A 5' FISH Probe covers some genomic sequences adjacent to the 5' (start) portion of the KMT2A gene. LSP KMT2A 3' FISH Probe covers the 3' (end) part as well as sequences downstream of the gene. LSP MLLT1 5' FISH Probe covers the 5' (start) portion of the MLLT1 gene and some adjacent genomic sequences. LSP MLLT1 3' FISH Probe covers the 3' (end) part as well as sequences downstream of the gene. The probe set is optimized to reveal translocations between the two genes.

Cat. No.	Volume
CT-PAC182-10-OG	10 Tests (100 μL)

Signal Pattern Interpretation

Normal Patterns **Abnormal Patterns** 202G* Other Patterns

*Overlapping orange and green signals can appear as yellow.

V2024.01.01

¹⁾ Tkachuk DC, et al. *Cell*. 71(4):691-700 (1992). 2) Huret JL, et al. *Leukemia*. 7(2):152-60 (1993). 3) Andreasson P, et al. *Eur J Haemato*l. 65(1):40-51 (2000). 4) Archimbaud E, et al. *Leukemia*. 12(1):25-33 (1998). 5) Barber KE, et al. *Genes Chromosomes Cancer*. 41(3):266-71 (2004).

CytoTest Inc. **IVD** 1395 Piccard Drive, Suite 308 Rockville, MD 20850, USA

^{*} CE IVD only available in certain countries. All other countries are either ASR or RUO. Please contact your local dealer or our headquarters for more information. DCN032