

ENGLISH

For Professional Use Only

DDIT3/CCP12 Tri-color Break Apart FISH Probe Kit

Introduction


The DDIT3/CCP12 Tri-color Break Apart FISH Probe Kit is designed to detect rearrangements in the human DDIT3 gene, located on chromosome band 12q13.3, along with the number of chromosome 12 copies per cell. In addition to revealing breaks, which can lead to translocation of parts of the gene, inversion, or its fusion to other genes, the probe set can also be used to identify other DDIT3 aberrations such as deletions or amplifications. Rearrangements and abnormal expression of the DDIT3 gene – also known as CHOP, CEBPZ, CHOP10, CHOP-10 or GADD153 - have been observed in myxoid liposarcoma and other malignancies.

Intended Use

To detect rearrangements in the human *DDIT3* gene located on chromosome band 12g13.3.

Cont.	Color
LSP DDIT3 5' FISH Probe	CytoGreen
LSP DDIT3 3' FISH Probe	CytoOrange
CCP12 FISH Probe	CytoAqua

Probe Design

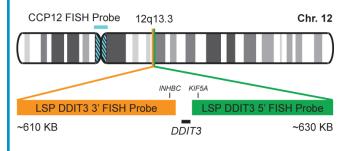

LSP DDIT3 5' FISH Probe covers some genomic sequences adjacent to the 5' end of the DDIT3 gene. LSP DDIT3 3' FISH Probe covers some sequence downstream of the 3' end of the gene. The two probes are flanking sequences across the DDIT3 gene in which variable breakpoints have been observed. CCP12 FISH Probe, derived from chromosome 12-specific alpha satellite DNA, is designed to serve as a control to determine the number of chromosome 12 copies per cell.

Not to Scale

Cat. No.	Volume
CT-PAC058-10-GOA	10 Tests (100 μL)

Signal Pattern Interpretation

Normal Patterns **Abnormal Patterns** 202G2A* Other Patterns

*Overlapping orange and green signals can appear as yellow.

¹⁾ Panagopoulos I, et al. *Oncogene*. 12(3):489-94 (1996). 2) Panagopoulos I, et al. *Biochem Biophys Res Commun*. 279(3):838-45 (2000). 3) Pérez-Losada J, et al. *Oncogene*. 19(52):6015-22 (2000). 4) Pérez-Mancera PA, et al. *Oncogene*. 21(11):1679-84 (2002). 5) Rabbitts TH, et al. *Nat Genet*. 4(2):175-80 (1993).

CytoTest Inc. **IVD** 1395 Piccard Drive, Suite 308 Rockville, MD 20850, USA

^{*} CE IVD only available in certain countries. All other countries are either ASR or RUO. Please contact your local dealer or our headquarters for more information.